## Uputstvo za upotrebu

VOX vaga PW 2001

Tehnoteka je online destinacija za upoređivanje cena i karakteristika bele tehnike, potrošačke elektronike i IT uređaja kod trgovinskih lanaca i internet prodavnica u Srbiji. Naša stranica vam omogućava da istražite najnovije informacije, detaljne karakteristike i konkurentne cene proizvoda.

Posetite nas i uživajte u ekskluzivnom iskustvu pametne kupovine klikom na link:

<https://tehnoteka.rs/p/vox-vaga-pw-2001-akcija-cena/>

- GBR
- ITA
- DEU
- RUS
- SRB
- ESP
- PRT
- GRC
- SVN
- HRV
- BIH
- MNE
- MKD
- ALB

**PW 2001**

OPERATING INSTRUCTIONS BODY-ANALYZER SCALE  
MANUALE DI ISTRUZIONI BILANCIA PER L'ANALISI CORPOREA  
BEDIENUNGSANLEITUNG KORPERANALYSEWAAGE  
ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ ВЕСЫ-АНАЛИЗАТОР ТЕЛА  
UPUTSTVO ZA UPOTREBU VAGA ZA ANALIZU TELA  
MANUAL DE USUARIO BÁSCULA DE ANÁLISIS CORPORAL  
MANUAL DO USUÁRIO BALANÇA DE ANÁLISE CORPORAL  
ΟΔΗΓΙΕΣ ΧΡΗΣΗ ΖΥΓΑΡΙΑ  
NAVODILA ZA UPORABU TEHTNICA ZA ANALIZO TELESA  
UPUTE ZA UPORABU VAGA ZA ANALIZU TIJELA  
UPUTSTVO ZA UPOTREBU VAGA ZA ANALIZU TIJELA  
УПАТСТВА ЗА РАКУВАЊЕ ВАГА ЗА АНАЛИЗА НА ТЕЛОТО  
MANUALI I PERDORUESIT PESHORJA PER ANALIZEN E TRUPIT

[www.voxelectronics.com](http://www.voxelectronics.com)

# Body-Analyzer Scale Instructions

PW 2001

Thank you for choosing our digital Bluetooth scale. The scale offers data of your personal weight, body fat percentage, body water percentage, muscle mass, bone mass, BMI (Body Mass Index) and more. The scale Provides you an accurate measurements of your overall health and fitness.

To ensure continuous operation and accuracy, please read these operating instructions carefully before use. By passing certain procedures may result in incorrect readings. We hope that you enjoy using your new Bluetooth scale.

### What's included

Your Bluetooth Scale box includes:

- Bluetooth smart scale x 1pc
- AAA batteries x 3pc
- Instruction manual x 1pc

### IMPORTANT SAFETY PRECAUTIONS:

When using digital scale, certain basic safety precautions must be followed, including the following:

1. Scale is a high precision measuring device. Never jump or stomp on the scale or disassemble. Please handle the scale carefully to avoid breakage.

1

2.To ensure weighing accuracy, when the scale is first used, please place them on an even surface and wait at least 10 seconds before stepping on the scale. It also applies if the scale has been moved. We suggest you to taking your measurements at the same time every day to ensure the accuracy.

3.Check the device before each use. Do not use the device or an electrode if it is damaged in any way. The continuous use of a damaged unit may cause injury, improper results, or serious danger.

4.Please check the batteries if the scale malfunction. Change to new batteries if needed.

5.When it is not in use, please remove the batteries and place the device in a dry room and protect it against extreme moisture, heat, lint, dust and direct sunlight. Never place any heavy objects on the equipment.

6.The scale contains plastic materials, if you are Allergic to plastic, do not use it.

For this Bluetooth scale, there are other safe precautions must be followed, including the following:

- 1.Use the Scale while barefoot. The scale cannot measure body fat percentage if you're wearing shoes or socks.
- 2.Make sure your feet are dry before weighing yourself.
- 3.Place the scale on a hard, flat surface during use. Using the scale on a soft surface - like carpet or linoleum - may affect the accuracy of the Scale readings.

2

4.The Bluetooth scale should not be used to diagnose or treat any medical condition. You should consult a doctor when undertaking any diet or exercise program. If you are pregnant/pacemaker or other internal device holder, you should not rely the data provided by this product. Please note that all the data are for reference and is not a substitute for advice from a licensed health-care professional.

5.This scale is not waterproof; never immerse the scale in water. To clean the surface, use a damp cloth or a glass cleaner. Do not use soap or other chemicals.

### FEATURES & SPECIFICATIONS:

- 1.This scale use BIA (Bioelectrical Impedance Analysis) technology to measure body composition data by sending a safe low-tension electrical current through the body.
- 2.High precision strain gauge sensors system
- 3.Auto ON/OFF
- 4.Weight Units: kg, lb (switch the units by following OPERATING INSTRUCCION 2.0)
- 5.Weight limit: 5kg-180kg(11lb-396lb)
- 6.Graduation increments: 0.05 kg/0.1lb
- 7.Measures Body Weight, B.M.I.(Body Mass Index), Body fat, Muscle, Water, Vessel Fat( V-fat), Bone Mass, Metabolism, Protein, Obesity Degree, Body age, Weight without fat.

3

- 8.Male/Female option
  - 9.Maximum 8 users number (add user members by following OPERATING INSTRUCCION 1.9)
  - 10.Overload indication: "EEr". Low batteries indication: "Lo" When the battery runs out, the display will show "Lo". Please replace new batteries.
  - 11.Power supply: AAA batteries.
- The typical service life of the new and unused batteries for AAA batteries is 4months measurements for the operation time for AAA batteries is 480 times.


### OPERATING INSTRUCTIONS:

1. General
  - 1.1 Remove all packaging before use.
  - 1.2 Remove insulating strip beneath the batteries in the battery holder at the back on the scale.
  - 1.3 Place the scale on hard, flat ground.
  - 1.4 Download the app: Bitafit from App store /Google Play and open mobile Bluetooth.


4

1.5 Register your own ID by email for keeping record. Add personal data during register and confirm OK.

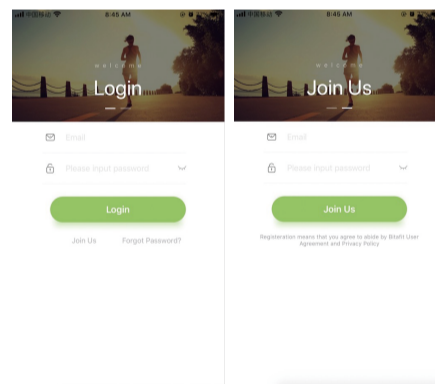
1.6 Turn on Bluetooth system on your phone, press 'Go' to add the scale.

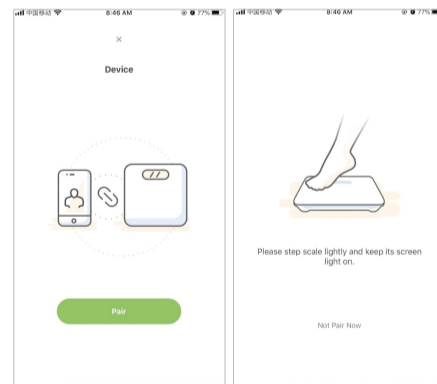
5

1.7 Step on the scale with BARE feet and wait for 10-15 seconds, the circle is stable with "ding" sound when you turn on bell ring, it means weighing and analyzing are finished.

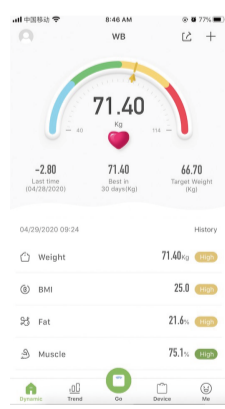
1.8 Click each data, you can find the analyzing report on Weight,BMI, Fat, Muscle, Water, V-fat, Bone Mass, Metabolism age. When you click the "chart" item at bottom, you can see your long-term tracking in chart form.

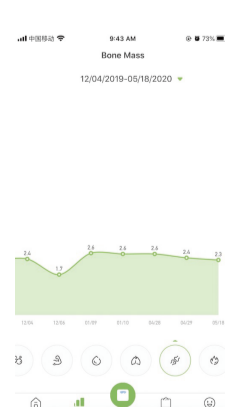
6

1.9 When you enter my account, you can add unlimited user member by press "add" and also add new users' personal data.

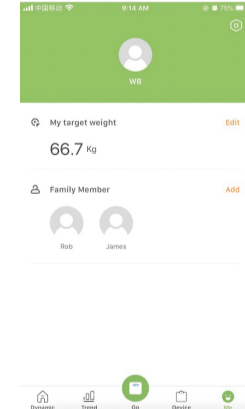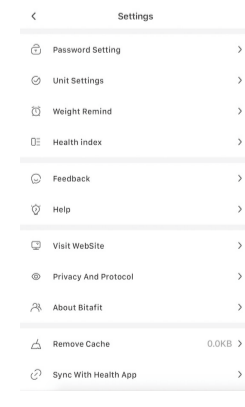2.0 Themes, Set goal, Weight Unit, Weight done sound, Password, and etc can be set on My Account page.

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